


# DB Orthodontics

## Winter Offers

Offers valid until February 28<sup>th</sup>

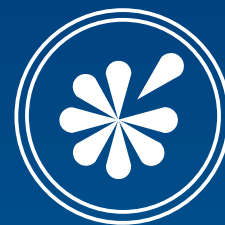

### Mini LED Curing Light

**NEW!**

High power density of 3,000 W/cm<sup>2</sup> when using 5.5mm booster tip. Takes only 4 seconds to bond stainless steel orthodontic brackets, and 2 - 3 seconds to bond ceramic brackets.

- LCD screen displays time selection
- Louder audible bleep or 4 second flash
- Compact and lightweight
- Silent - No Fan
- Light guide rotates 360°
- Standby indicator and low battery warning

Cordless - Noiseless - Lightweight

DB06-0065  
NRP £870.00

**Special Offer £750.00**

7 Day Trial  
Available



Freephone 0800 783 3552



mini implant system

**Sterile Implants Now Available**



Colour Coded	Diameter (mm)	Body Length	Neck Length	Sterile Code	Typical Insertion Sites
	1.5	6	Short	<b>DB10-0009</b>	Anterior Maxilla - perpendicular insertions Mandible - perpendicular insertions
	1.5	9	Short	<b>DB10-0010</b>	Anterior Maxilla - buccal alveolus Posterior Maxilla - all insertions Mandible - oblique insertions
	1.5	9	Long	<b>DB10-0012</b>	Maxilla - palatal alveolus
	2.0	6	Long	<b>DB10-0015</b>	Mid-palate
	2.0	9	Long	<b>DB10-0016</b>	Edentulous areas Temporary restorative abutment

NRP £43.95  
**Launch Offer £35.15**

**Sterile Complete Kits**



- Contains;
- Infinitas sterilisation tray
  - Standard handle
  - Standard inserts; Screwdriver, Mucotome, Bone Punch
  - 5 mini implants of your choice
  - Mini handle
  - Mini inserts: Screwdriver, Mucotome, Bone punch

NRP £875.00  
**Launch Offer £700.00**

**Infinitas Courses**



**Course Dates**

- London 26th Feb 2010
- Leeds 12th March 2010

**Objectives of the course;**

- To critically appraise the literature evidence on mini-implant techniques.
- To explain the design, benefits and clinical usage of the Infinitas mini-implant system.
- To provide practical in vitro teaching on the use of the Infinitas system.

**New Advanced Course**

- Leeds 11th March 2010

Further enquires contact DB Orthodontics on 0800 783 3552 or

**Book Online At: [infinitas-miniimplant.co.uk](http://infinitas-miniimplant.co.uk)**



The evolve stainless steel bracket system features a two piece bracket, metal injection moulded from advanced 17-4 low nickel alloy brazed to an 80 gauge mesh base.

The bracket slot is CNC milled and finished with an advanced polishing process resulting in a perfectly smooth slot.

The mushroom shaped hook secures chains and attachments and prevents slipping in between patient appointments.

Also available with gingivally offset bicuspids brackets.

Available in .022 Roth and MBT.

**From 60p Per Bracket**



Quantity Rates (NRP)	Price / Kit	Price /Bracket
1 case kit	<b>£30.00</b>	<b>£1.50</b>
10 case kits	<b>£25.00</b>	<b>£1.25</b>
50 case kits	<b>£19.80</b>	<b>£0.99</b>
100 case kits	<b>£12.00</b>	<b>£0.60</b>

**Ultra Low Friction Slot**

Advanced polishing process creates mirror finish by removing high surface peaks at molecular level.

**Beveled Tie Wings**

To reduce visibility from the front without compromising on strength.

**Rounded Tie Wings**


For Patient Comfort

**Centre Scribe Line**


For accurate placement in line with the axis of the tooth

**Palladium Alloy Base**

Assures strongest possible brazed joint and resists electrolytic reaction that can cause brazing failure during use.

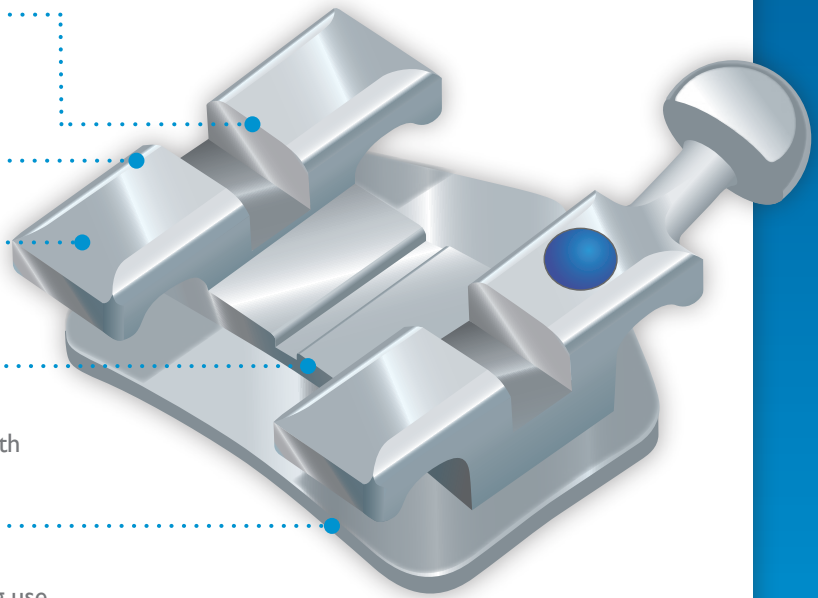

ID marks with enhanced colours for easy identification.

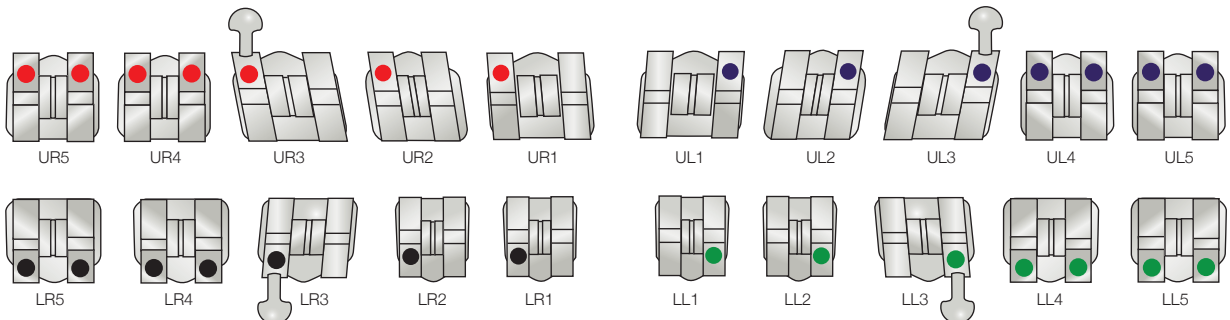

# Buy 2 Get 1 Free\*



College Tweezers  
DB05-0319  
**£7.55**



Reverse Keat Tube Holder  
DB05-0423  
**£21.65**



Flat Plastic  
DB05-0316  
**£8.30**



Probe  
DB05-0312  
**£11.30**



Mitchell Trimmer No. 4  
DB05-0317  
**£10.55**



Divider  
DB05-0  
**£13.40**



6" Ruler  
DB05-0318  
**£8.20**



Mosquito Forceps - Standard  
DB05-0315  
**£8.25**



Mathieu Needle Holder T.C. Tips  
DB05-0406  
**£36.25**



Mershon Band Pusher  
DB05-0559  
**£9.95**

## Mirrors No. 4

Tested to withstand 2000 cycles in an autoclave. The design of the reflecting layers provides maximum image quality and extremely high resistance to deterioration. Heads packed in boxes of 12. Handles packed individually.

Mirror heads	DB05-0313M -	<b>£13.15/pk12</b>
Mirror handle	DB05-0313 -	<b>£1.75 each</b>
Mirror & handle	DB05-0313H -	<b>£2.85 each</b>



\* Cheapest item free